





VentecLife.com

To place orders:

Email: ventec-orders@reacthealth.com

Telephone: 1-844-698-6276

Ventec Life Systems, a wholly owned subsidiary of React Health, Inc., reserves the right to make changes in specifications, pricing, and/or to discontinue any product at any time without notice or obligation and will not be liable for any consequences resulting from the use of this publication.

Caution: U.S. Federal law restricts these devices to sale by or on the order of a physician.

MKT-00004 Rev L
Copyright © 2024 Ventec Life Systems



Configurable Software

In addition to the multiple hardware configurations, VOCSN software can be configured to disable and enable available therapies based on patient needs.

Multi-Function Ventilator

CUSTOMIZABLE TO EACH PATIENT

VOCSN is one model with multiple configurations to meet your patient needs and care environments. Get all five therapies or just the mix of therapies needed.

Configuration	Ventilation	Internal O2 Concentrator	Cough	Suction	Nebulizer	High Flow Nasal Cannula Therapy	High Pressure External O2 and FiO2 Monitor	Low Pressure External Oxygen
V*Home (PRT-01445-000)	✓	—	—	—	—	✓	—	✓
V+Pro (PRT-01198-000)	✓	—	—	—	—	✓	✓	✓
VC (PRT-01125-000)	✓	—	✓	—	—	✓	—	✓
VC+Pro (PRT-00973-000)	✓	—	✓	—	—	✓	✓	✓
VCSN+Pro (PRT-01172-000)	✓	—	✓	✓	✓	✓	✓	✓
VOCSN (PRT-01126-000)	✓	✓	✓	✓	✓	✓	—	✓
VOCSN+Pro (PRT-00853-000)	✓	✓	✓	✓	✓	✓	✓	✓



V*Home



V+Pro



VC+Pro



VOCSN

CIRCUIT TYPE & OPTIONS

Switch therapies without changing circuits

The Ventec One-Circuit allows patients to switch between therapies with the touch of a button. Rather than several cumbersome tubes connecting to separate devices, the Ventec One-Circuit eliminates clutter to create one easy-to-manage system for better care, comfort, and peace of mind.

Active

The streamlined design and dual bellows of the active circuit provide the performance of an active valve with the size and convenience of a passive valve. The patent pending valve design supports cough and ventilation functions in a single circuit, and also combines the drive and sense lines into a single trilumen tube. Additionally, VOCSN uses advanced software technology to help prevent historical problems associated with condensation to optimize Vte sensing, accuracy, and triggering by purging condensation on every breath.

Passive

The Ventec passive circuit is light, small, and easy to use. The unique passive valve is designed to close the leak during cough exsufflation, thus making Touch Button Cough possible without changing circuits.

Valveless

The valveless Ventec One-Circuit is compatible with masks that include a leak and high flow therapy.

Oxygen Direct®



The Oxygen Direct system minimizes oxygen lost to leaks and dead space, improves battery life, and makes using supplemental oxygen easier. VOCSN includes advanced pulsing technology and an independent O2 tube inside the patient circuit. The system delivers oxygen directly to the patient interface during the first two thirds of the patient's breath, when the patient needs oxygen most. Note: For use with VOCSN only.

Heated

The Ventec One-Circuit heated in-wall wire system provides consistent heat distribution throughout the patient circuit to ensure comfortable humidification.

Exhalation Filter Cuff

The exhalation filter cuff covers the exhalation valve on the Ventec One-Circuit, and provides a connection point for a bacterial filter to mitigate the risk of cross-contamination from exhaled pathogens.

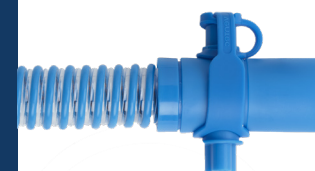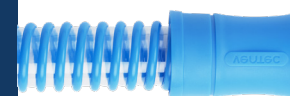

ACTIVE



Active

Filter cuff
on adult circuits

Active valve

Elastic circuit tube clips

6 ft. length

10-PACK

Adult (w/cuff) PRT-00950-001
Pediatric PRT-01011-000



Active, Heated

Filter cuff
on adult circuits

Active valve

In-wall heated wire

Elastic circuit tube clips

5 ft. length

10-PACK

Adult (w/cuff) PRT-00947-001
Pediatric PRT-01012-000



Active, Oxygen

Filter cuff
on adult circuits

Active valve

Integrated Oxygen Direct tube


Elastic circuit tube clips

6 ft. length

Note: For use with VOCSN only

10-PACK

Adult (w/cuff) PRT-00930-001
Pediatric PRT-01013-000

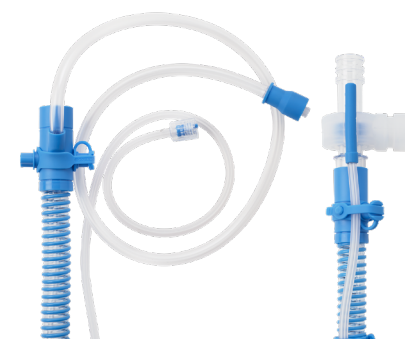

Active, Heated, Oxygen

Filter cuff
on adult circuits

Active valve

In-wall heated wire

Integrated Oxygen Direct tube and
extension for use with humidifier

Elastic circuit tube clips

5 ft. length

Note: For use with VOCSN only

10-PACK

Adult (w/cuff) PRT-00945-003
Pediatric PRT-01014-001

PASSIVE



Passive

Passive valve

6 ft. length



Passive, Heated

Passive valve

In-wall heated wire

Elastic circuit tube clips

5 ft. length



Passive, Oxygen

Passive valve

Integrated Oxygen Direct tube

6 ft. length

Note: For use with VOCSN only



Passive, Heated, Oxygen

Passive valve

In-wall heated wire

Integrated Oxygen Direct tube and extension for use with humidifier

Elastic circuit tube clips

5 ft. length

Note: For use with VOCSN only

10-PACK

Adult	PRT-00946-000
Pediatric	PRT-01015-000

10-PACK

Adult	PRT-00948-000
Pediatric	PRT-01016-000

10-PACK

Adult	PRT-00949-000
Pediatric	PRT-01017-000

10-PACK

Adult	PRT-00944-001
Pediatric	PRT-01018-001

VALVELESS



Valveless

Valveless circuit for use with vented mask

6 ft. length



Valveless, Heated

Valveless circuit for use with vented mask. Pediatric size is intended to be used for high flow therapy for pediatric and adult patients

Elastic circuit tube clips

5 ft. length


10-PACK

Adult	PRT-00942-000
Pediatric	PRT-00983-000

10-PACK

Adult	PRT-00943-000
Pediatric	PRT-00985-000

	Active	Passive	Valveless
Non Heated	10-PACK Adult (w/cuff) PRT-00950-001 Pediatric PRT-01011-000	10-PACK Adult PRT-00946-000 Pediatric PRT-01015-000	10-PACK Adult PRT-00942-000 Pediatric PRT-00983-000
Heated	10-PACK Adult (w/cuff) PRT-00947-001 Pediatric PRT-01012-000	10-PACK Adult PRT-00948-000 Pediatric PRT-01016-000	10-PACK Adult PRT-00943-000 Pediatric PRT-00985-000
Oxygen	10-PACK Adult (w/cuff) PRT-00930-001 Pediatric PRT-01013-000	10-PACK Adult PRT-00949-000 Pediatric PRT-01017-000	
Heated, Oxygen	10-PACK Adult PRT-00945-003 Pediatric PRT-01014-001	10-PACK Adult (w/cuff) PRT-00944-001 Pediatric PRT-01018-001	

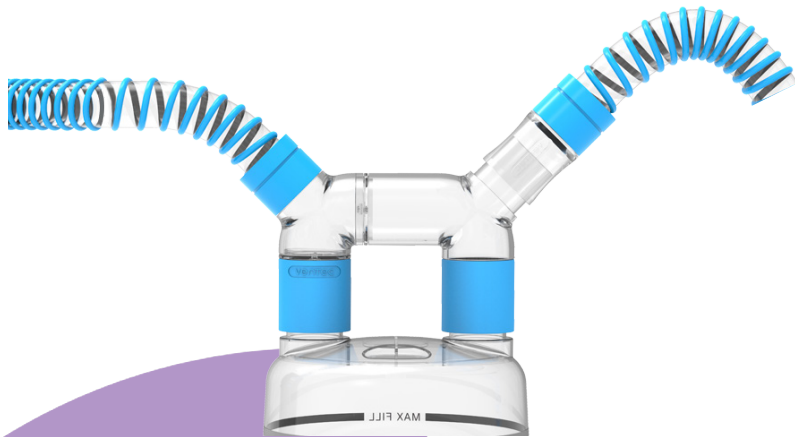

Humidifier Bypass

The humidifier bypass allows the patient to remain connected to the humidifier while using the Touch Button Cough therapy with a heated Ventec One-Circuit. The bypass attaches to the humidifier water chamber to allow humidified air to reach the patient during ventilation and insufflation. Then, during exsufflation, check valves redirect the flow of air to bypass the water chamber of the humidifier (humidifier chamber not included)


10-PACK
PRT-00989-000



O2 Low-Pressure Inlet Adapter

The O2 low-pressure inlet adapter inserts into the low-pressure oxygen port on the back of VOCSN to connect to a low-pressure oxygen source.

10-PACK
PRT-00981-000

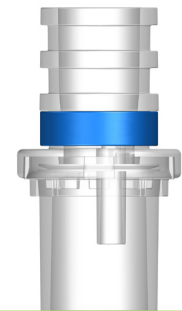

Passive Valve

The passive valve provides a fixed leak exhalation port during ventilation, with a 1-way valve designed to close during the exsufflation phase of cough assist therapy

10-PACK
PRT-01451-000



External Bacterial Filter

10-PACK

PRT-00850-000

The primary, low resistance, external bacterial filter captures 99.99% of bacteria and viruses. Easily attaches directly to VOCSN and the Ventec One-Circuit. More than 50% smaller in size than traditional filters to make transport easier and reduce the risk of disconnects



Internal Bacterial Filter

10-PACK

PRT-00976-000

The secondary, low resistance, internal bacterial filter captures 99.99% of bacteria and viruses to protect against cross-contamination. Serviceable by trained technicians, typically during 10,000 hour scheduled maintenance



Fan Filter

25-PACK

PRT-00978-000

Washable and reusable foam filter cleans air entering the cooling fan intake



Patient Air Intake Filter

25-PACK

PRT-00977-000

Washable and reusable foam filter cleans air entering the patient air intake



Travel Suction Canister

300 mL detachable travel suction canister attaches to VOCSN to provide suction therapy on the go. Includes an internal filter that will self seal to protect VOCSN if the travel suction canister overfills

10-PACK
PRT-00849-000



External Suction Canister Adapter

External suction canister adapter connects VOCSN to any third party external suction canister. Includes an internal filter that will self seal to protect VOCSN if the external suction canister overfills

36 in. length


5-PACK
PRT-01002-000



Suction Canister Rollstand Holder

Suction canister rollstand holder easily mounts to the VOCSN Rollstand center pole or drop pole accessory to securely hold third party suction canisters up to 1,200 mL. Convenient clip holds suction tubing (suction canister and tubing not included)

1-PACK
PRT-01010-000

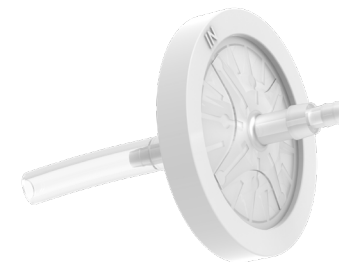

Nebulizer Filter

The VOCSN nebulizer filter is an antibacterial filter for use with nebulizer therapy to minimize cross-contamination between uses, with a BFE filtration efficiency of >99.9%. The filter will also self-seal to prevent fluid in the nebulizer cup from reaching VOCSN

5-PACK
PRT-01092-000



Everyday mobility

VOCSN provides everyday mobility for patients with up to nine hours of on board battery (ventilation only) with the internal battery and two removable, rechargeable batteries.



VOCSN AC Power Adapter

Two-prong AC power cable and AC/DC power adapter. The power connector includes a threaded locking mechanism to reduce the risk of accidental disconnects. Wire retention clip ensures the power supply cord remains firmly connected to the adapter. LED on the power adapter indicates proper power connection. The power supply also comes with two hook-and-loop straps to manage excess cabling. This power adapter recharges the three lithium-ion batteries while VOCSN is on or off.

PRT-01333-000



24V Wheelchair Power Cable

Use to power VOCSN with 24V wheelchair outlets. One end of the cable includes the VOCSN power connector with threaded locking mechanism to reduce the risk of accidental disconnects. The charger output connector on the other end of the 72" cable uses a 3-pin XLR male plug. This power adapter recharges the three lithium-ion batteries while VOCSN is on or off.

PRT-01027-000



Removable, Rechargeable Battery Set, Lithium-ion

Set of two lithium-ion 14.4V 5.8Ahr hot-swappable batteries with integrated push button charge indicator. Batteries can be charged in VOCSN or with the detachable battery quick charger (battery quick charger sold separately).

PRT-00980-000



Removable Battery Quick Charger

60W desktop charger for two VOCSN removable, rechargeable batteries. Indicator lights display battery charge status. Each battery charges fully in three hours. The charger is powered by a standard two-prong AC cord (batteries not included).

PRT-01003-000

VOCSN Roll Stand

The VOCSN roll stand is rugged, mobile, and designed for easy accessorizing. Each roll stand is assembled in two easy steps and includes:



PRT-00569-001



Nurse Call Cable: Open (Black)

Connects VOCSN to normally open remote alarm, nurse call, or other alarm systems that sense contact closure through a 1/4 in. (.6 cm) phono jack; 10 ft. coiled cable

10-PACK PRT-00851-000



Nurse Call Cable: Closed (White)

Connects VOCSN to normally closed remote alarm, nurse call, or other alarm systems that sense contact closure through a 1/4 in. (.6 cm) phono jack; 10 ft. coiled cable

10-PACK PRT-00852-000



Carry Bag

The VOCSN carry bag supports everyday mobility with an included shoulder strap and multiple attachment points to secure VOCSN to a wheelchair, bed rail, or other mount. The Carry Bag is designed to accommodate full VOCSN functionality while attached. It is compatible with the roll stand, travel suction canister, and allows access to both removable, rechargeable batteries

PRT-01143-000



React DataLink™

The React DataLink provides wireless support for patient data sharing, allowing clinicians to view trends, track patient compliance with prescribed therapy, and facilitate informed treatment decisions. Access to connected device data is managed on a dedicated, password protected web portal called Multi-View Connect™, access to which requires cellular network subscription.

PRT-01446-000

USB Drive

USB 2.0 drive, with FAT32 format, compatible with all VOCSN configurations

5-PACK PRT-01416-000



Recommended Accessories

These accessories are available from third-party manufacturers and are compatible for use with VOCSN.



Humidifiers

The Fisher & Paykel HC550 (Home Care) and MR850 (Hospital) humidifiers are compatible with the Ventec humidifier bypass

VentecLife.com/humidifier550
VentecLife.com/humidifier850



6 L/min Nebulizer Cup

The Sunset Disposable Nebulizer Kit with T-Piece is recommended for use with nebulizer therapy to ensure accurate VOCSN flow compensation



Wheelchair Carrier

The Freedom Wheelchair Carrier was custom-designed to attach VOCSN to a wheelchair

VentecLife.com/carrier



Mouthpiece Circuit Kit

The Mouthpiece Circuit Kit connects to VOCSN to provide patients with mouthpiece ventilation

VentecLife.com/mouthpiece



← Seattle & Univ. of Washing

→ Airport & SeaTac/Airport →

← Seattle & Univ. of Washington → Airport & SeaTac/Airport →



REACTHEALTH